

Saflo 90 Catheter




The Saflo 90 is a complementary product to the original Saflo and provides a 90 degree catheter configuration for full control of depth of insertion. It has an effective passive needle safety system in compliance with the sharps safety regulations.

The principal features of this set are:

- Soft catheter infusion set for maximum patient comfort and safety
- Catheter at right angles to tape fixation for consistent infusion depth
- Translucent cap for infusion site visibility
- Passive safety system completely protects the user from the dangers of sharps injury
- Available in 6 mm and 9 mm catheter lengths for optimum patient compatibility
- Colour coded for ease of identification of needle length
- Easy insertion and removal of catheter
- Low priming volume connector tubing (less than 0.2mL)
- Suitable for a wide range of clinical applications
- Soft catheter cap for maximum patient comfort

Saflo 90 Port

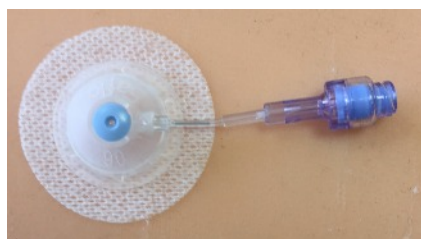

With the same basic characteristics as the Saflo 90 range the new port configuration adds extra flexibility. It may be used as a means of administering intermittent therapies or, in conjunction with standard administration sets, for higher flow applications such as subcutaneous rehydration therapy.

Saflo 90 Order Codes

Code	Box Qty	Description	Cap
SF90-106	12	Saflo 90 6mm Catheter 100cm Tubing	Green
SF90-109	12	Saflo 90 9mm Catheter 100cm tubing	Blue
SF90-109-AS	12	Saflo 90 9mm Catheter 100cm tubing with anti-siphon valve	Blue
SF90-P09	12	Saflo 90 9mm Needle free infusion port	Blue
SF90-P06	12	Saflo 90 6mm Needle free infusion port	Green