


Saflo 90 Catheter- Immunology Applications



The principal features of this set are:

- Soft catheter infusion set for maximum patient comfort and safety
- Catheter at right angles to tape fixation for consistent infusion depth
- Translucent cap for infusion site visibility
- Passive safety system completely protects the user from the dangers of sharps injury
- Available in 6 mm and 9 mm catheter lengths for optimum patient compatibility
- Colour coded for ease of identification of needle length
- Easy insertion and removal of catheter
- Low priming volume connector tubing (less than 0.2mL)
- Suitable for a wide range of clinical applications
- Soft catheter cap for maximum patient comfort

Saflo 90 Port

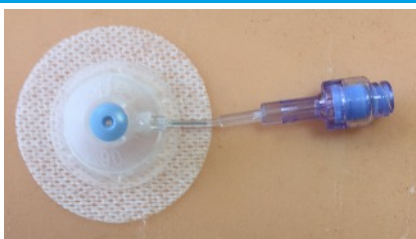

With the same basic characteristics as the Saflo 90 range the port configuration adds extra flexibility. It may be used as a means of administering intermittent therapies or, in conjunction with standard administration sets, for higher flow applications such as subcutaneous rehydration therapy.

The construction provides a fluid path equivalent to that of a 24g needle which is suitable for use with more viscous solutions such as Immunoglobulins. This may be coupled with an appropriate extension set to provide a flexible connection system for a variety of applications.

Currently the Saflo range offers a choice of a single or two way (Y- Configuration) extension set for use in conjunction with the Saflo 90 port to provide this flexibility. A 3-way configuration is currently under development.

