

Saflo Catheter



The Saflo Catheter is the heart of the Saflo subcutaneous infusion system. It is a completely new design and is the first subcutaneous device specifically produced to meet the requirements of all current sharps safety legislation.

By using a passive safety system the Saflo meets all the requirements of this legislation and of the relevant British Standard.

The system is latex free and has a PVC free fluid path to ensure the widest compatibility with common drugs. Low priming volume tubing is used throughout.

Connector Tubing



This provides the basic unit for connection between the pump or infusion source and the catheter inserted into the patient's tissue.

It incorporates a simple, removable, latching system for connection to the catheter and is available in tubing lengths of 100 cm and 60 cm. Also available with anti-siphon valve.


Connector Tubing with Side Port



The side port is provided for cases where it is appropriate to inject additional medication alongside the standard infusion, for example, in the treatment of breakthrough pain.

The needle-free injection port provides the safest means of administration for the extra medication and eliminates the need for additional needle or catheter insertion into the patient.


Injection Port

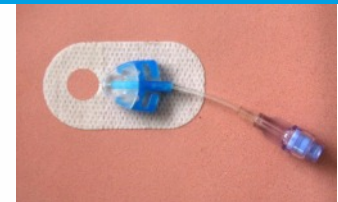

This option is provided for cases where it is only necessary to provide occasional injections or prior to use of the catheter for longer term therapies. The system has a priming volume of approximately 0.1 mL and is not considered suitable for insulin administration.

It may be appropriate to flush the infusion line after drug administration, in accordance with local hospital policy.

Saflo 90 Catheter

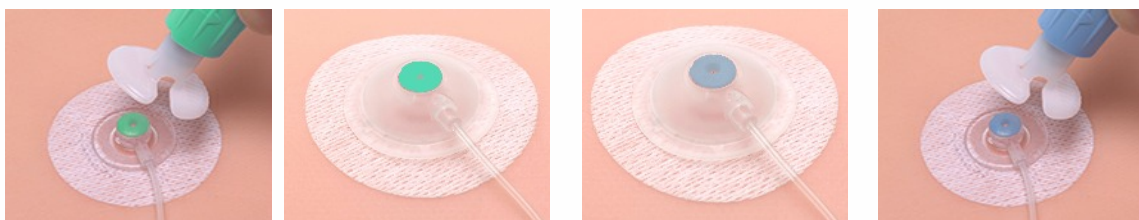

The Saflo 90 is a complementary product to the original Saflo and provides a 90 degree catheter configuration for full control of depth of insertion. It has an effective passive needle safety system in compliance with the sharps safety regulations.

The principal features of this set are:

- Soft catheter infusion set for maximum patient comfort and safety
- Catheter at right angles to tape fixation for consistent infusion depth
- Translucent cap for infusion site visibility
- Passive safety system completely protects the user from the dangers of sharps injury
- Available in 6 mm and 9 mm catheter lengths for optimum patient compatibility
- Colour coded for ease of identification of needle length
- Easy insertion and removal of catheter
- Low priming volume connector tubing (less than 0.2 mL)
- Suitable for a wide range of clinical applications
- Soft catheter cap for maximum patient comfort
- Some codes available with integral anti-siphon valve

Saflo 90 Port



With the same basic characteristics as the Saflo 90 range the port configuration adds extra flexibility. It may be used as a means of administering intermittent therapies or, in conjunction with standard administration sets, for higher flow applications such as subcutaneous re-hydration therapy.

Available in both 6 mm and 9 mm catheter lengths.

Saflo Microvolume Extension sets



The Saflo Extension sets use the same tubing construction as the general Saflo range. This is polyethylene lined PVC tubing with a priming volume of approximately 0.2 mL.

This may be used with any of the Saflo products to extend the line length for extra flexibility.

They are available in 100 cm and 150 cm lengths with optional anti-siphon

Saflo System Order Codes

Order Code	Description
SFC001-00	Catheter Introducer Assembly
SFT100-00	Saflo connector tubing 100 cm
SFT100-SP	Saflo connector tubing 100 cm with side port
SFP001-00	Saflo infusion port connector
Microvolume Extension Set Codes	
SFX100-AS	Microvolume extension tubing 100 cm with anti-siphon valve
SFX150-00	Microvolume extension tubing 150 cm
SFX-102Y	Microvolume Extension set Y- configuration

All above SF codes are packed in boxes of 10 units

Saflo 90 Order Codes

Order Code	Description	Cap Colour
SF90-106	Saflo 90 6 mm Catheter 100 cm Tubing	Green
SF90-109	Saflo 90 9 mm Catheter 100 cm tubing	Blue
SF90-P06	Saflo 90 6 mm Needle free infusion port	Green
SF90-P09	Saflo 90 9 mm Needle free infusion port	Blue

All SF90 codes are packed in boxes of 12 units